



My C.L.E.

CINETIC

LIFT

EXPERT

**DISCOVER THE POWER OF ADVANCED ENERGIES
THE 1ST AT HOME BEAUTY DEVICE***

LEDs AND MICRO-CURRENT TECHNOLOGIES
TO SHAPE THE FIRMNESS, TEXTURE, CLARITY
AND VITALITY OF THE SKIN.




TAKE YOUR BEAUTY ROUTINE
TO THE NEXT LEVEL :

My C.L.E.

With My C.L.E., inspired by Carita's professional beauty machine, the Cinetic Lift™ Expert, upgrade your daily beauty routine, at home.

- Featuring two types of energy to stimulate your skin, day after day - Micro-currents* and LED** lights intensely energize the skin from within to improve its texture, density and vitality.
- Multiple uses and benefits provide a bespoke beauty routine - My C.L.E. can be individually adjusted to give your skin exactly what it needs at any given time.

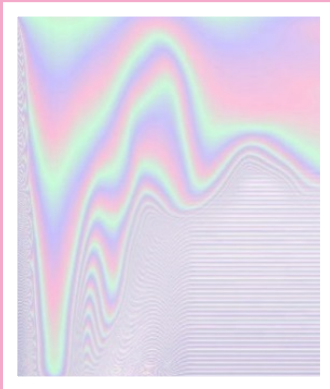

MICRO-CURRENTS*:

« **THE INSTANT LIFTING ENERGY** »
AVAILABLE IN 5 LEVELS OF
INTENSITY, IT ENHANCES VOLUME
AND TONE OF SKIN



LEDs**:

« **THE PROGRESSIVE TAILOR-MADE ENERGY** » AVAILABLE IN 4 COLORS FOR
AN ACTION ON FIRMNESS - TEXTURE -
CLARITY - VITALITY



FIRMNESS



TEXTURE



CLARITY



VITALITY

VISIBLE AMPLIFIED RESULTS

Tested and approved by women.

**VISIBLE RESULTS
FROM 4 WEEKS**

FACE BEAUTY PROGRAMS *

MORE THAN **80%**

Feel that their skin after 28 days is
TRANSFORMED, REDENSIFIED, AS LIFTED

79%

Feel that their skin after 28 days is
MORE RESISTANT, YOUNGER

83%

Are convinced by the **EFFICACY** of the routine and would like to add My C.L.E. to their daily routine

BODY BEAUTY PROGRAMS**

96%

Feel that their skin after 28 days is
MORE COMFORTABLE

71%

Feel that their skin after 28 days is
SMOOTHER

71%

Feel that the texture of the skin after 28 days is
MORE HOMOGENEOUS

*PROTOCOL: RED LED AND LEVEL 5 MICRO-CURRENTS INTENSITY WITH LE SERUM NEOMORPHOSE 7MN/DAY DURING 28 DAYS.

**SELF-ASSESSMENTS ON 49 WOMEN WHO APPLIED THE CARITA CRÈME RÉVÉLATRICE FERMÉTÉ IN ASSOCIATION WITH MY C.L.E. ONCE A DAY : ON THIGHS, STOMACH AND ARMS (RED LED LIGHT AND MICRO-CURRENTS INTENSITY 5) - RESULTS AFTER 4 WEEKS

AT HOME FACE BEAUTY PROGRAMS:

ANTI-AGEING: READY, SET, LIFT!



RED LED & MICRO-CURRENTS INTENSITY 4-5

3 SERUM OPTIONS

FOR A FILLING AND ANTI-WRINKLE ACTION

LE SÉRUM NÉOMORPHOSE

(6-8 pressures for full face)

FOR A TIGHTENING AND FIRMING ACTION

LE SÉRUM STIMULIFT

(6-8 pressures for full face)

FOR A GLOBAL ANTI-AGEING ACTION


LE SÉRUM PARFAIT 3 ORS

(4-6 pressures for full face)



MASSAGE YOUR FACE WITH UPWARD LIFTING MOTIONS TO FIRM AND PLUMP YOUR SKIN. REPEAT 3 TIMES EACH MOVEMENT. START THE PROGRAM FROM THE CHEST UP TO THE FOREHEAD.

PURENESS: FIGHT YOUR IMPERFECTIONS!

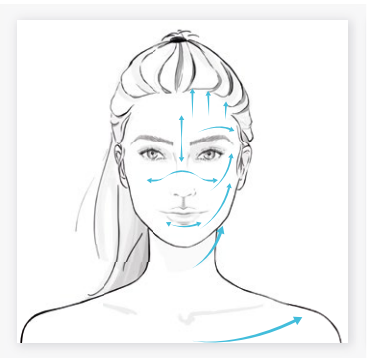

BLUE LED & MICRO-CURRENTS INTENSITY 1-3

1 SERUM OPTION

FOR A BALANCED AND PURIFIED SKIN ACTION

LE SÉRUM CONTRÔLE


(6-8 pressures for full face)



MASSAGE YOUR FACE WITH GENTLE BACK AND FORTH MOVEMENTS LOCATED ON T-ZONE (CHIN, NOSE, FOREHEAD) TO PURIFY AND REBALANCE YOUR SKIN. REPEAT 3 TIMES EACH MOVEMENT. START THE PROGRAM FROM THE CHEST UP TO THE FOREHEAD.

MAKE SURE TO ALWAYS CLOSE YOUR EYES WHEN YOU USE THE DEVICE ON THE CHEEKS, TEMPLES AND FOREHEAD

ANTI-REDNESS/DARK SPOT: MISSION FLAWLESS TONE!

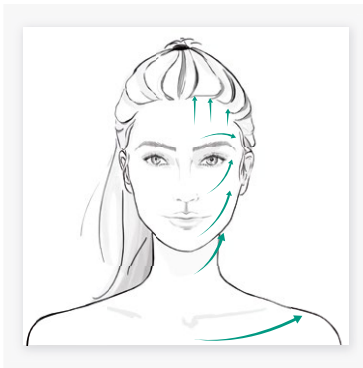

GREEN LED & MICRO-CURRENTS INTENSITY 1-3

3 SERUM OPTIONS

FOR AN EVEN TONE AND ANTI-DARK SPOTS ACTION

LE SÉRUM CLARTÉ
(6-8 pressures for full face)

FOR A SOOTHING AND ANTI-DISCOMFORT ACTION

LE SÉRUM SENSIDOTE
(6-8 pressures for full face)

FOR A GLOBAL ANTI-AGEING ACTION


LE SÉRUM PARFAIT 3 ORS
(4-6 pressures for full face)



MESSAGE TARGETED AREAS USING CIRCULAR MOTIONS TO PURIFY YOUR SKIN AND REVEAL A MORE EVEN AND HOMOGENOUS COMPLEXION.

REPEAT 3 TIMES EACH MOVEMENT. START THE PROGRAM FROM THE CHEST UP TO THE FOREHEAD.

RADIANCE: SUPERGLOW ON THE GO!

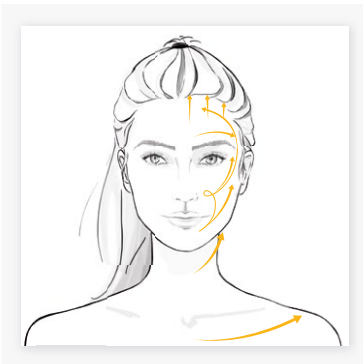

WHITE LED & MICRO-CURRENTS INTENSITY 1-3

2 SERUM OPTIONS

FOR A REHYDRATING AND ANTI-DULLNESS ACTION

LE SÉRUM LAGON
(6-8 pressures for full face)

FOR A GLOBAL ANTI-AGEING ACTION

LE SÉRUM PARFAIT 3 ORS
(4-6 pressures for full face)



MESSAGE YOUR FACE WITH A «3» SHAPED MOVEMENT TO REVITALIZE YOUR SKIN AND REVEAL A RADIANT COMPLEXION.

REPEAT 3 TIMES EACH MOVEMENT. START THE PROGRAM FROM THE CHEST UP TO THE FOREHEAD.

OUR 5 COMMITMENTS

1 - TAILOR RESULTS

My C.L.E. HAS BEEN TESTED ON ALL SKIN TYPES AND OFFERS A TAILOR-MADE SOLUTION. FIND THE PROGRAM THAT BEST SUITS YOUR SKIN'S NEEDS, WHETHER YOU WANT TO TARGET FIRMNESS, SKIN TEXTURE, EVENNESS OR RADIANCE.

2 - TESTED AND APPROVED

My C.L.E. IS RECOMMENDED BY ALL OF OUR BEAUTY EXPERTS.

3 - EASY TO USE


FIND COMPREHENSIVE HELP AND TUTORIAL VIDEOS ON OUR WEBSITE (WWW.CARITA.FR) WITHIN My C.L.E. SECTION.

4 - EASY TO CLEAN

ALL YOU NEED IS A COTTON BALL AND A BIT OF WATER OR ALCOHOL. My C.L.E. MUST NOT HAVE DIRECT CONTACT WITH WATER.

5 - 2-YEAR WARRANTY

FOR FURTHER INFORMATION PLEASE CONSULT THE WEBSITE (WWW.CARITA.FR).

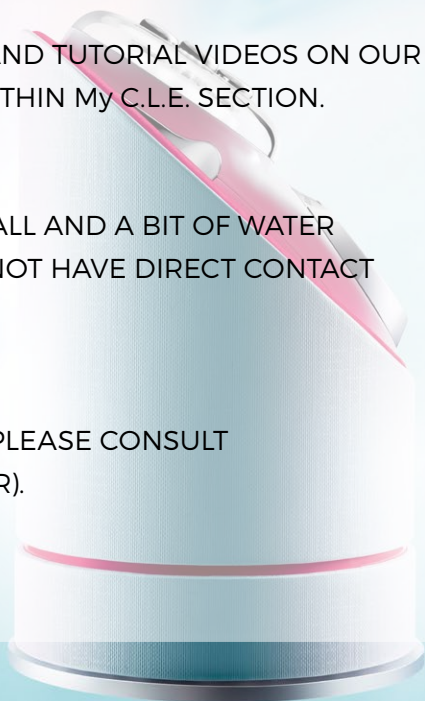




CARITA
PARIS




AT HOME BODY BEAUTY PROGRAMS:

TOP ARMS CHALLENGE




AFTER A MANUAL MASSAGE WITH LA CRÈME RÉVÉLATRICE FERMETÉ ON THE ARMS, PERFORM DEEP SMOOTHING MOVEMENTS FOLLOWED BY A STATIC MOVEMENT FOR 3 SECONDS. FINISH YOUR MASSAGE WITH PINCHING CIRCULAR MOVEMENTS FOR A MORE TARGETED ACTION. THE SKIN IS SMOOTHED AND TONED.


1. APPLY THE FORMULA

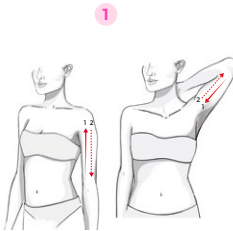

X3

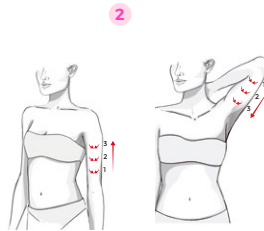


X3

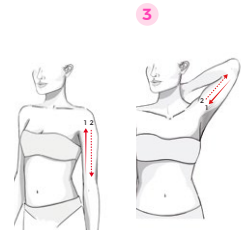


X3


2. ACTIVATE

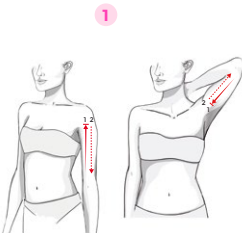

X3

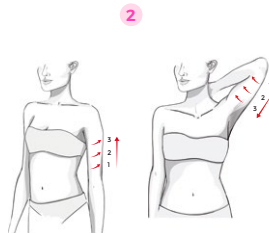

X3

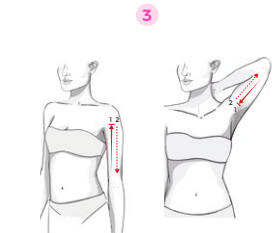

X3

SMOOTHING MOVEMENTS
PINCHING CIRCULAR MOVEMENTS

LA CRÈME RÉVÉLATRICE FERMETÉ
(1 TO 2 PEA-SIZED DABS)

RED LED & MICRO-CURRENTS INTENSITY 5

WAIST:
SLOW CIRCULAR MOVEMENTS
PINCHING CIRCULAR MOVEMENTS

ARMS & LEGS:
STATIC LIFT MOVEMENTS
PINCHING CIRCULAR MOVEMENTS

MAKE SURE TO ALWAYS KEEP THE DEVICE A HAND'S WIDTH AWAY FROM THE SENSITIVE AREA AROUND THE ARMPITS

OBJECTIVE SUBLIME WAIST


5
MIN

AFTER A MANUAL MASSAGE WITH LA CRÈME RÉVÉLATRICE FERMETÉ ON THE STOMACH, PERFORM SLOW CIRCULAR MOVEMENTS FOLLOWED BY PINCHING CIRCULAR MOVEMENTS. THE SKIN IS SMOOTHED AND GRADUALLY TIGHTENED.

1. APPLY THE FORMULA



X3

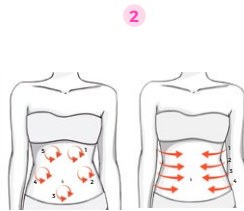


X3




X3


2. ACTIVATE

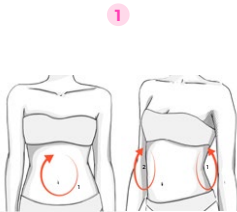

X3

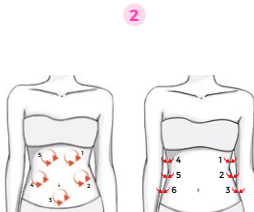

X3

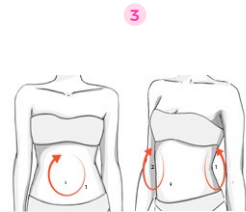


X3

MISSION FLAWLESS LEGS



AFTER A MANUAL MASSAGE WITH LA CRÈME RÉVÉLATRICE FERMETÉ ON THE THIGHS, PERFORM DEEP SMOOTHING MOVEMENTS FOLLOWED BY A STATIC MOVEMENT FOR 3 SECONDS. FINISH YOUR MASSAGE WITH PINCHING CIRCULAR MOVEMENTS FOR A MORE TARGETED ACTION. DAY AFTER DAY, THE SKIN LOOKS FIRMER.

1. APPLY THE FORMULA

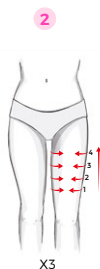

2. ACTIVATE

